

VISION SCREENER



VISIOCLICK®

AUTOMATION AT YOUR FINGERTIPS

AUTOMATION OF SCREENING FOR VISUAL DISORDERS

- ▶ Adapts to your Visiolite® equipment
- ▶ Automation of all tests
- ▶ Multilingual vocal instructions
- ▶ Compact & lightweight





PERFORMANCE IN THE LEAD

Convert your Visiolite® Master or Master GT to an automated device. VisioClick also adapts to previous versions (from 2012). Compact, lightweight, and **operational from the first connection**. A simple update of the Visiolite® software is all it takes. No operator presence required during tests. Multilingual voice instruction ensures global understanding.

STEP 1



Patient is equipped with the headphones and receives vocal instructions for tests displayed on the screen.

STEP 2



Patient clicks on response unit to answer instructions.

STEP 3




Results are automatically saved.

STEP 4



Information is given to the patient and operator indicating the test is completed. Results are recorded and can be interfaced with all «business software».

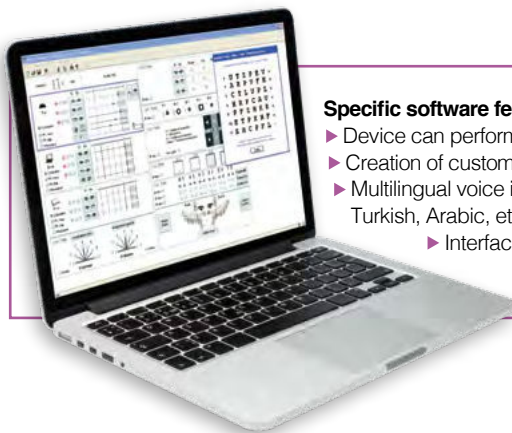

VISIOLITE® SOFTWARE

Specific software features:

- ▶ Device can perform all tests individually in the order of your choice or a pre-programmed entire test protocol
- ▶ Creation of customizable profession-related sequences that can be performed in automatic mode
- ▶ Multilingual voice instructions (French, English, German, Spanish, Dutch, Italian, Polish, Portuguese, Turkish, Arabic, etc.)
- ▶ Interfacing with market software packages (MEDTRA, HEALTH@WORK, STETHO, PREVENTIEL, KITRY, DINAMIT, PADOA...)

YOUR DISTRIBUTOR



**1 FREE BOX
OF 50 PAIRS OF
EARPAD COVERS***

Single-use biocompatible
for all your tests.

Earpad covers are recommended.

Benefits:

- Increases lifespan of earpads
- A more hygienic solution
- 11 cm diameter fits a wide range of headphones
- No sound reduction

* VisioClick® delivered with a pack of earpad covers.

ACCESSORIES

(Included with your device)

- The Headphones + stand
- 1 travel case
- 1 box of biocompatible earpad covers

2
YEARS
WARRANTY
VISIOCLICK®

1
YEAR
WARRANTY
HEADPHONES



51 rue Antoine Primat - CS 60194
69625 Villeurbanne Cedex France

Tél +33 (0)4 72 34 89 89
Fax +33 (0)4 72 33 43 51
contact@fim-medical.com

